## INVACARE® MATRX®

# PSP (Posture Seat Polymer)



The **Invacare Matrx PSP** (Posture Seat Polymer) cushion combines the stability and support of the PS foam base with a polymer overlay for extra protection and shear reduction. An ultra-lightweight cushion providing zero maintenance, worry free skin protection and positioning.

**Call to order** 

**1.800.333.6900** invacare.com

**1.888.433.6818** motionconcepts.com








HCPCS Code: Pediatric/Standard - E2607 Heavy-Duty - Pending (E2608)

#### **Features and Functional Benefits**

- Startex<sup>™</sup> outer cover Reversible can be used smooth or fabric side out
- Moisture resistant and breathable to keep cool and dry
- Includes hook and loop strips
- Startex<sup>™</sup> inner cover for extra moisture protection
- Molded high resilient polyurethane foam made with Ultra-fresh™ anti-microbial
- Rear cushion radius and waterfall front edge
- Gently contoured to provide neutral leg position
- Molded polymer insert for extra protection and shear reduction
- Unique waffled relief zone allows immersion in ischial area
- Designed to provide off-loading to hips and thighs for support and stability

#### **Optional Accessories**

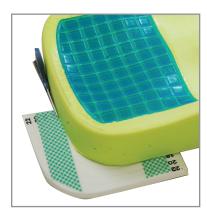

**POK** – Pelvic Obliquity Kit (2 pieces) accommodates up to one inch of obliquity



**Cushion Rigidizer** – lightweight and low profile















Model #	Width	Depth	Height	Height at Leg Abductor/Adduct	Product tor Trough	Weight	Limited Capacity
Warranty							
PSP (Pediatric)	10" to 14"	10", 12", 14", 16"	3.5"/3.5"	2.75"	1.2 lb. and up	150 lb.	2 Years
PSP (Standard)	16", 17", 18", 20"	16", 18", 20"	3.5"/3.5"	2.75"	2.2 lb. and up	300 lb.	2 Years
PSP (Heavy-Duty)	22", 24", 26"	18", 20", 22"	3.5"/3.5"	3.00"	4.6 lb. and up	600 lb.	2 Years







### Invacare Corporation

www.invacare.com

#### USA

One Invacare Way Elyria, Ohio 44035-4190 (800) 333-6900

#### Canada

570 Matheson Blvd. E. Unit 8 Mississauga, Ontario L4Z 4G4 Canada (800) 668-5324

#### **Motion Concepts**

 $www.motionconcepts.com\\info@motionconcepts.com$ 

#### USA

700 Ensminger Rd., Suite 112 Tonawanda, NY 14150 (888) 433-6818 Fax: (888) 433-6834

#### Canada

84 Citation Dr., Unit 1-7 Concord, ON L4K 3C1 (866) 748-7943 (905)-695-0134 Fax: 905-695-0138

#### © 2021

Invacare Corporation. All rights reserved. Trademarks are identified by the symbols <sup>™</sup> and ®. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Matrx is a registered trademark of Motion Concepts. Otto Bock is a trademark of Otto Bock HealthCare. Specifications are subject to change without notification. Form No. 18-296 rev. 1/12 160959