

vive health



VIVE BARIATRIC ROLLATOR MOB1040


The Vive bariatric rollator walker provides stabilizing support while walking and a comfortable seat when resting. This rollator is made with a durable steel frame with height-adjustable handles, a padded backrest, and an 18" wide padded seat for greater comfort. It can also support up to 500 pounds. The bariatric rollator folds flat for easy portability and includes an under-the-seat bin for storing everyday items.

WHAT'S INCLUDED

- Vive Bariatric Rollator
- Removable storage bin

COLORS  

EASILY MANEUVER INDOORS AND OUT

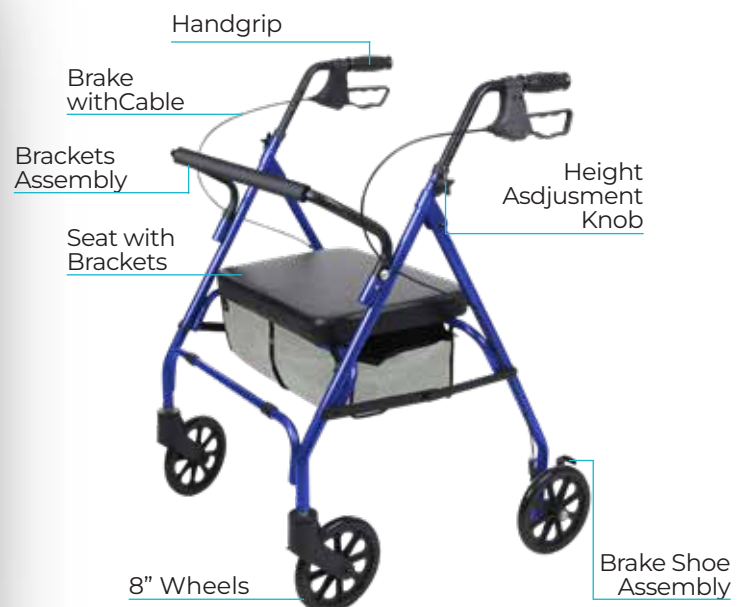

www.aracent.com



(888) 550-7750



sales@aracent.com



Specifications

Base when open 28.5" by 20" by 38"

Width between the handles 22"

Adjustable Height 34" to 38"

Seat 18" by 13.5"

Seat height from the floor 22"

Wheels 8"

Storage bin 20" x 10" x 5"

Materials Steel frame

Storage Bag Dimensions 15" x 13" x 3.5"

Storage Bag Weight Capacity 11 lbs

STABLE MOBILITY SUPPORT

Providing stabilizing support, the Vive bariatric rollator allows you to easily navigate indoors and out and comfortably rest when needed. The heavy-duty rollator is wider than standard rollators for a comfortable walking experience and features a comfortably wide seat and padded backrest with an under-the-seat storage bin for safely carrying essential items.

STRONG STEEL FRAME

Safely supporting up to 500 pounds, the Vive bariatric rollator is made with a durable steel frame that requires minimal, tool-free assembly. The heavy-duty rollator handles are also height-adjustable between 34" to 38" for a customized fit.

FOLDS FOR EASY STORAGE

Designed for easy transportation, the bariatric rollator walker folds flat for convenient storage. The durable rollator weighs approximately 27 pounds.

LOCKING BRAKING SYSTEM

Each of the ergonomic handles is angled to position hands in a neutral position, reducing fatigue and allowing easy access to the bicycle-style hand brakes. The loop brakes also allow you to lock the wheels in place, ensuring the rollator will not move when activated.



SR22- 9R1



www.aracent.com



(888) 550-7750



sales@aracent.com